




JT-17 Arch Trimmer Operation Manual

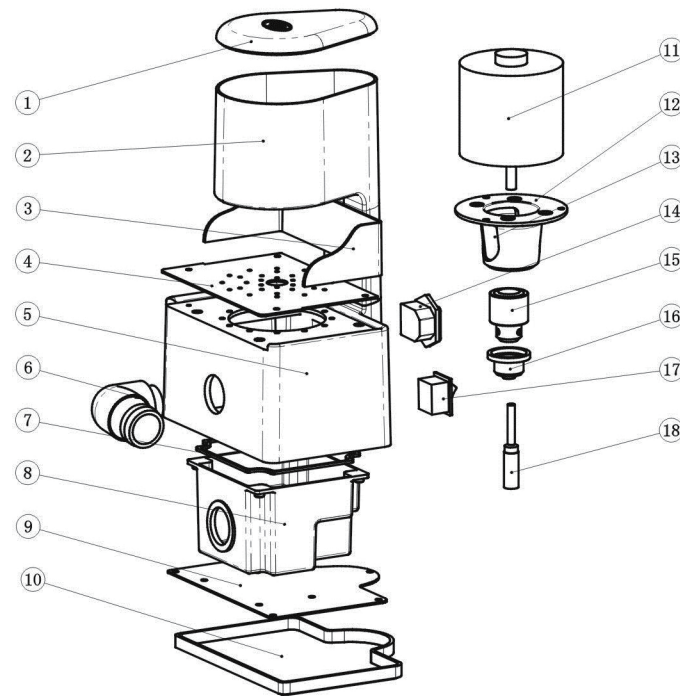

No.	Part List	Q' TY
1	upper cover board	1
2	upper part	1
3	dust-proof plastic frame	1
4	workboard	1
5	lower part	1
6	dust-collecting bending tube	1
7	sealing rubber ring	1
8	built-in powder-collector box	1
9	bottom plate	1
10	rubber bottom plate	1
11	motor	1
12	bearing sleeve	1
13	screw cap	1
14	socket	1
15	chuck	1
16	Rubber sleeve	1
17	power switch	1
18	rotatable file	1

Power supply :220v/50hz (110v/60hz)

Power :100Watt

Rotate Speed: 2800~3000rpm

1. Please check whether it is in normal operation after switch on the power. If not, please check whether the built-in fuse of the socket are fusing. If yes, please promptly replace it (10A);
2. the front-end of the rotatable file must be below the work table 0.5 ~1mm, ensuring effective grinding;
3. If you want to replace the rotatable file or adjust its height, please turn the dust mask by 90 °, with two attachment wrenches hold the upper spindle and below respectively and lock the outside edge of the screw cap tightly to make it tight or loose. After work, recover the dust mask again.
4. Set the dust-collecting tube into the mouth of the dust rowing to work together with dust-collector, and dust section with cleaner. If long-term work cause excessive dust within the box, Pls discharge the bottom plate firstly, then the powder box and clean it.
5. please pay attention to the safety of operation.
 - 1). never touch the rotatable file when it is on work;
 - 2) never operate it with gloves.

Contact Us

E-mail Address:

Dentalsalemall.com@gmail.com

Telephone Number: +86-023-68438829

Official Site:<https://www.dentalsalemall.com>