
Dentium Advanced Sinus Kit

DASK | Sinus Elevation Drill, Stopper & Hand Instrument Set

KIT LAYOUT

SCHEMATIC VIEW - NOT TO SCALE

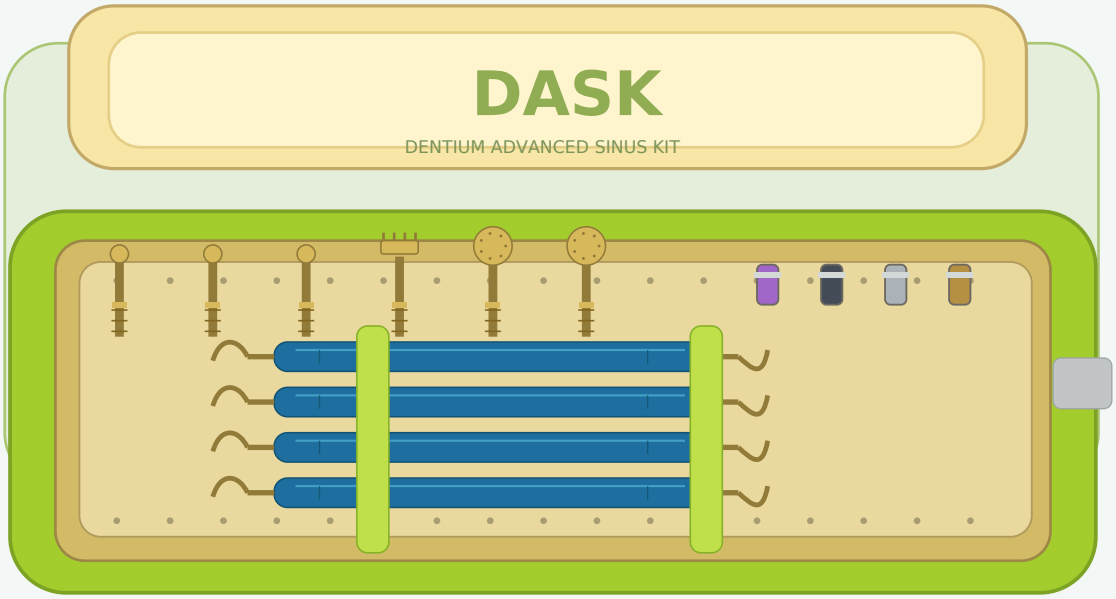

06 DRILLS

04 DRILL STOPS

04 HAND INSTRUMENTS

PRODUCT OVERVIEW

The DASK instrument set is designed for dental sinus elevation procedures using crestal and lateral approaches. The set combines approach-specific drills, drill stops and sinus elevation instruments in an organized tray.

01 CRESTAL APPROACH

Three dedicated drills for crestal access.

02 LATERAL APPROACH

Three dedicated drills for lateral access.

03 DRILL STOPPERS

Four depth-specific stopper options.

PRODUCT IDENTIFICATION

Model / Reference: DASK

Country of origin: Korea (as listed by the seller)

Components & Specifications

DRILL SET / 6 PIECES

REFERENCE NUMBER & IRRIGATION

DRILL	APPROACH	PART NO.	IRRIGATION
#1	Crestal	XRT332035	Internal
#2	Crestal	XRT372035	Internal
#3	Crestal	XED331035D	Internal
#4	Lateral	XRT064025	Internal
#5	Lateral	XRT084025	Internal
#6	Lateral	XST083025	External

DRILL STOPPERS / 4 PIECES

FOR DRILLS #1-3

STOPPER	DRILL DEPTH	LENGTH (L)	PART NO.
08	8 mm	10.6 mm	XFDST08
06	6 mm	12.6 mm	XFDST06
04	4 mm	14.6 mm	XFDST04
02	2 mm	16.6 mm	XFDST02

SINUS ELEVATION HAND INSTRUMENTS / 4 PIECES

XSE1L

XSE2L

XSE3L

XSE4L

RECOMMENDED DRILLING PARAMETERS

DASK Drills #1-6	**800-1,200 rpm**	**30-45 Ncm**

Irrigation: internal on drills #1-5; external on drill #6.

USE & CONFIGURATION NOTES

For trained dental professionals. Surgical indications, instrument selection, drilling settings, cleaning and sterilization must follow the applicable Dentium instructions for use. Confirm supplied components and reference numbers before clinical use.

Contact Us
E-mail Address:
Dentalsalemall.com@gmail.com
Telephone Number: +86-023-68438829
Official Site:https://www.dentalsalemall.com