Algimax-Ⅱ

GX300


Alginate Mixer



Smooth mixture in seconds for a quick and accurate impression

Patented anti-vibration mechanism

Maximize workflow and hygiene

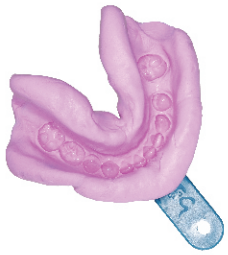

Algimax-Ⅱ GX300 Alginate Mixer


The ALGIMAX-II is an automatic alginate mixing machine designed for use in dental clinics and laboratories. It features a motor-driven mixing system with adjustable mixing times to help prepare alginate impression materials. The machine provides an automated mixing process that can help reduce manual mixing effort and support consistent material preparation. Mixing results may vary depending on the alginate material, powder-to-water ratio, and operating conditions.

Simple steps to creating a impression:



1. Add alginate powder and water into the mixing cup.

2. Moisten with the stirrer for a few seconds and close the mixing cup.

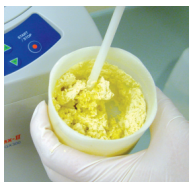




3. Put the mixing cup into the mixing well.

4. Cover the lid and press" START".



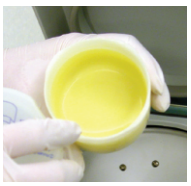


5. After mixing, use the spatula and remove the mixed impression.

3 Default Mixing Modes

Scoops	Mixing Time
	8 secs
	10 secs
	12 secs

Specifications:

- Timer Range 1 ~ 16 seconds
- Mixing Modes 3 (Default: 8, 10, 12 secs)
- Speed 3,400 rpm
- Power Input AC 110 or 220V, 50/60 Hz
- Dimensions 250(L) x 205(W) x 285(H) mm
- Net Weight 17.5 kg
- Accessories Spatula/Stirrer/Storage Case/ Water Dosing Bottle/Mixing Cups(x4)

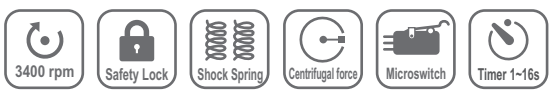

E-mail Address:
Dentalsalemall.com@gmail.com
Telephone Number: +86-023-68438829
Official Site:https://www.dentalsalemall.com

CE